

# King Vision® aBlade™

## Video laryngoscope



### Key Benefits

#### Affordable

- The King Vision® aBlade is an excellent solution to solve the issue of limited availability of advanced visualisation solutions in the OR.
- With King Vision®, you can improve patient care by having high performance visualisation capabilities for every intubation at an affordable price.

#### Portable

- King Vision® is lightweight, self-contained, battery operated and water resistant. King Vision® can be used continuously for about 80 minutes making it highly portable.
- The aBlades are individually packaged. Each blade is disposable, eliminating concerns of cross contamination. King Vision® offers the choice of two blade types depending on preferred tube delivery technique.

#### Durable

- King Vision® and aBlades have been designed to be your primary tool for all intubations.
- Video output capability is compatible with external monitor and recording devices.
- The robust, full-colour, non-glare display can resist repeated cleaning. The aBlade video adapter is reusable, designed with materials that can withstand repeated use.

## King Vision® video laryngoscope

Visionary performance in the palm of your hand. When it comes to visualising the airway for routine or difficult intubations, you need a tool that you can depend on each and every time. King Vision® is designed for indirect visualisation of the vocal cords which minimise the need for tissue manipulation and improves the ability to see anterior airways. Such reliable performance is a necessity, not an option. King Vision® comes with the choice of two blade types that allow users to keep their preferred technique in terms of tube delivery, offering a unique flexibility.

#### Full range of blade sizes

The King Vision aBlade system offers a complete range of video adapters and disposable blades for a wide range of patients. This allows you to have one solution to cover your intubation cases from infants to adults, whether in the OR, ICU or ER. The King Vision aBlade solution can also be implemented in the prehospital setting for Emergency Services.

#### Portable

The King Vision® display is lightweight, self-contained, battery operated, and water-resistant. The reusable display and video adapter come packaged in a protective case. aBlades are individually packaged making King Vision® very portable.

#### Affordable

The disposable blades allow for economical use of the King Vision® for most of your intubations. Combined with the high performance visualisation capabilities this offers great value while also delivering excellent patient care. Take the guess-work out of your intubations. See the vocal cords and confirm placement of the endotracheal tube.

#### Durable

The King Vision® is designed to be your primary tool for intubations. The high-quality display comes with a 1-year warranty. The robust, full-colour, smudge-proof, non-glare display can resist repeated cleaning and normal use wear and tear.



**FULL-COLOUR TFT LCD DISPLAY**  
Anti-reflective

**PROTECTIVE GASKET**

**ERGONOMIC DESIGN**  
Integrated blade/handle design for complete control

**CHANNELED AND STANDARD DISPOSABLE BLADES**  
for guided or standard insertions

**GUIDING CHANNEL**  
for easy ETT placement and removal

**ANTI FOG TREATMENT**  
in every blade

**EXCELLENT OPTICS**

Age	Size 1	Size 2	Size 3
0			
1			
2			
3			
4			
5			
6			
7			
8			
9			
10			
11+			

#### Size recommendation:

For infants and children the small video adapter 1-2 and associated blades are indicated:

- aBlade size 1 (Infant patients <4 years)\*
- aBlade size 2+2C (Child patients 1-10 years)\*


For larger children and adults the large video adapter 3-4 and associated blades are indicated:

- aBlade size 3+3C (Children >5 years and Adults)\*

\*Size ranges are based on teeth-to-vallecula data for the given ages, tracheal tube sizing recommendations and actual experience.

A medical professional must evaluate on a patient-by-patient basis.

#### aBlade sizes

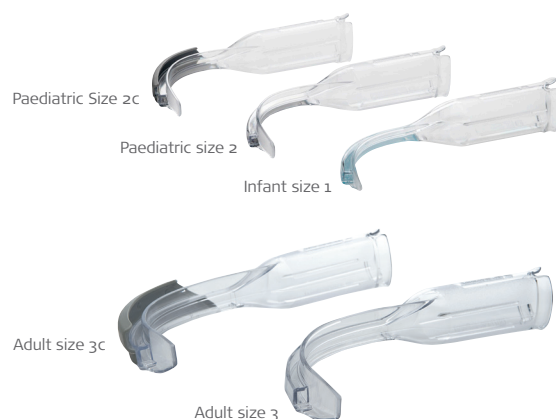

## Ordering information

Item no.	Description	Quantity
KVLAB1	Standard aBlade, size 1	20/CS
KVLAB2	Standard aBlade, size 2	20/CS
KVLAB2C	Channeled aBlade, size 2C	20/CS
KVLAB3	Standard aBlade, size 3	20/CS
KVLAB3C	Channeled aBlade, size 3C	20/CS
KVLVA12	aBlade Video Adapter 12	1/CS
KVLVA34	aBlade Video Adapter 34	1/CS
KVLABKIT3	aBlade Kit	1/CS
KVIS01VA34	1 Display and 1 Video Adapter 34 in aBlade case	1/CS
KVIS01ADCASE	aBlade case	1/CS
KVIS01	King Vision Display	1/CS

## Material information

Part	Material
Blade	Polycarbonate, TPE (Thermoplastic Polyurethane Elastomer)
Video adapter	Polycarbonate/ABS, TPE (Thermoplastic Polyurethane Elastomer)
Display	Polycarbonate/ABS anti-reflective coating on display window

### Sales Company

**Ambu Ltd.**  
8 Burrell Road  
St Ives  
Cambridgeshire  
PE27 3LE  
Tel 01480 498403  
Fax 01480 498405  
uksales@ambu.com  
www.ambu.co.uk

### Manufacturer

**King Systems**  
15011 Herriman Boulevard  
Noblesville, IN 46060  
Tel. 317 776 6823  
Tel. 800 642 5464  
Fax 317 776 6827  
www.Kingsystems.com

Distributed by  
**Ambu**

**CE** US: Rx only