



Key Features

- Transparent dome for visual check of bleeding, vomitus and spontaneous breathing
- Thumb rest for a tight, easy seal
- Comfortable design and texture
- The buttoned-on cuff enables easy removal for cleaning
- Can be autoclaved repeatedly at 134°C (except size 0A)
- Available in 7 sizes
- Size 0 with boring impedes ventilation pressures higher than approx. 40 cm H₂O, enabling the user to administrate correct ventilation on small children with a body weight of at least 15 kg (3 years)

Ambu® Silicone Face Mask

The Ambu Silicone Face Mask is designed for use with manual and automatic resuscitators and ventilators.

Ambu Silicone Face Mask is self-inflating for an easy and tight seal. The cuff pressure may be regulated (except for mask size 0A).

SIZE

0A Infant	28 g
0 Infant w. boring	39 g
0 Infant	38 g
2 Child	63 g
3/4 Adult/Child	78 g
4 Adult	80 g
5 Adult	105 g


MATERIALS

Cuff	Silicone rubber
Cuff size 0A	Silicone rubber with foam filling
Dome	Polysulphone
Mask connector	Silicone rubber



ORDERING INFORMATION

000 313 000	Ambu Silicone Face Mask size 0A Infant
000 314 000	Ambu Silicone Face Mask size 0 Infant w. boring
000 315 000	Ambu Silicone Face Mask size 0 Infant
000 316 000	Ambu Silicone Face Mask size 2 Child
000 312 000	Ambu Silicone Face Mask size 3/4 Adult/Child
000 317 000	Ambu Silicone Face Mask size 5 Adult
000 319 000	Ambu Silicone Face Mask size 6 Adult

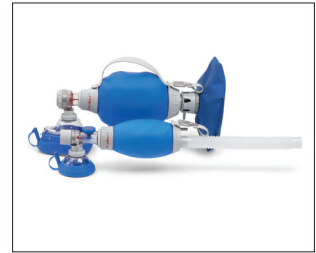

Ambu Silicone Face Mask connected to
Ambu Mark IV Resuscitator



Ambu Silicone Face Mask connected to
Ambu Silicone Resuscitator

Ambu A/S

Baltorpbakken 13
DK-2750 Ballerup
Denmark
T +45 72 25 20 00
F +45 72 25 20 50
www.ambu.com

