

## DATASHEET

---

Sterile straight from the pack, Ambu aScope Duodeno is a single-use endoscope with side viewing optics, deflectable tip and an elevator to control the position of compatible endoscopic accessories.

With aScope™ Duodeno, you have a new scope with consistent performance every time, and there is no risk of patient cross-contamination. What's more, there's no need for costly reprocessing or repair.



Distal tip showing elevator in use with a sphincterotome



Elevator control lever to raise (up) and lower (down) the elevator

# Ambu® aScope™ Duodeno Sterile. Single-use. Problem solved.



## KEY BENEFITS

---

Sterile

No risk of patient cross-contamination

Cost transparency

Eliminates the risk and cost of reprocessing

Familiar design and function to that of conventional duodenoscopes

# Ambu

## SPECIFICATIONS

### Optical System

Direction of View	Retro Side Viewing 6°
Field of View	130°

### Insertion Section

Distal End Outer Diameter	φ 13.7 mm
Insertion Tube Outer Diameter	φ 11.3 mm
Working Length	1240 mm

### Instrument Working Channel

Inner Diameter	φ 4.2 mm (12.6 FR)
----------------	--------------------

### Bending Section

Angulation Range	Up 120° / Down 90° Left 90° / Right 110°
------------------	---

### Illumination


Illumination source	Dual LEDs
---------------------	-----------

### Dimensions

Weight	700 g
--------	-------

## ORDERING SPECIFICATIONS

Reference number	Description	Quantity
485001000	aBox Duodeno	1
482001000	aScope Duodeno	Box of 6

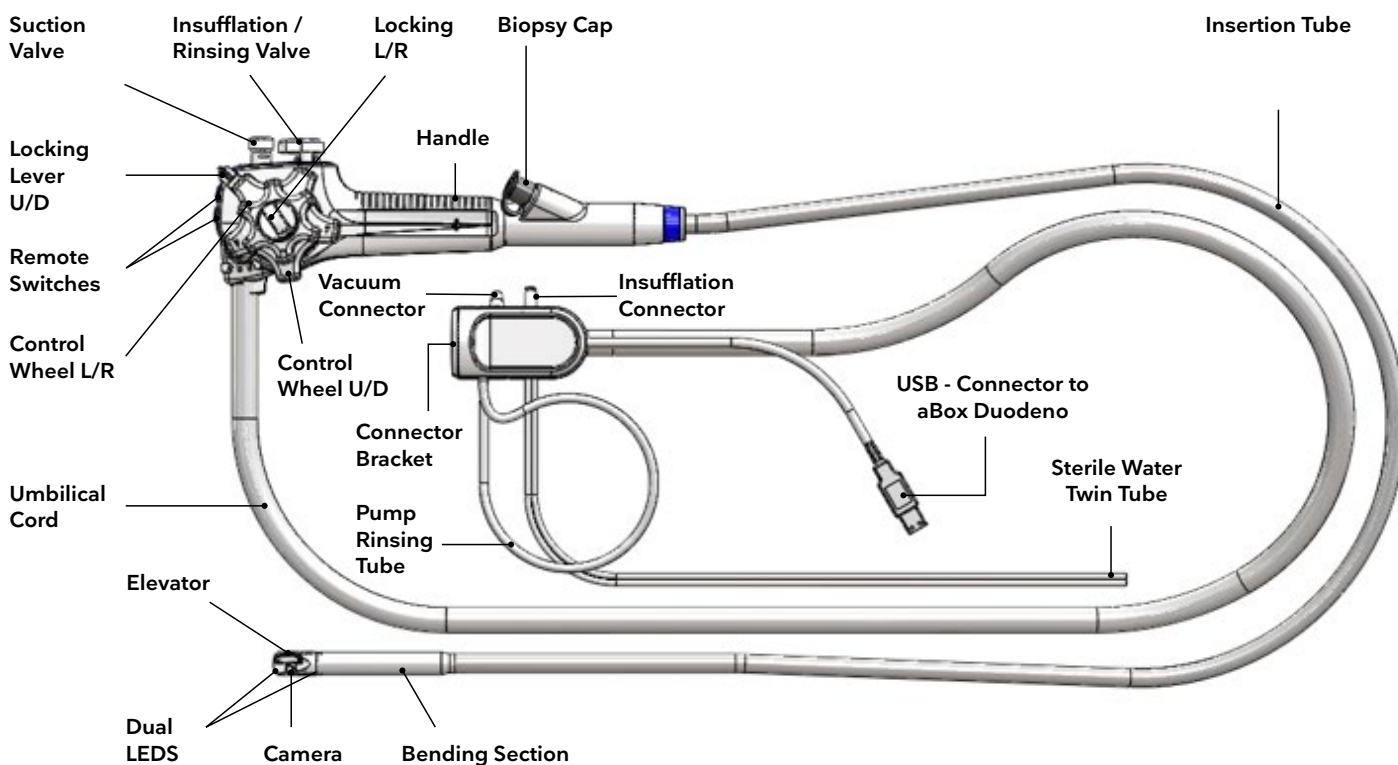

# Ambu

Ambu Ltd  
 Incubator 2  
 Alconbury Weald Enterprise Campus  
 Alconbury Weald  
 Cambridgeshire  
 PE28 4XA  
 uksales@ambu.com  
 www.ambu.co.uk



US: Rx Only